

## Warranty certificate Exidéal

This warrants the Product against defects in materials and workmanship under normal use for a period of one year from the date of retail purchase by the original distributor. Please fill out the form below the warranty with your personal information and purpose of usage.

Item's Stock Number		
Serial Number		
Warranty Term	<b>After purchasing for 1 year</b>	
Purchase Date		
Customer	ZIP	
	Address	
	Name	
	TEL	
Distributor		

### 【Regulation for free repair】

- 1) This product and its components will perform properly and reliably only when operated and maintained in accordance with the instructions contained in this manual and accompanying labels.
    - Please bring proof-of-purchase. Proof-Of-Purchase is required to verify warranty period at the original distributor. However, shipping charges for all warranty repairs are at the owner's expense when you send it to our customer service or original distributor.
    - Also business trip charges for all warranty repairs are at the owner's expense when outwork repair is requested.
    - Please contact our customer service if you decide not to ask for free-repair at the original distributor.
  - 2) Please contact our customer service about free-repair when you move to another place.
  - 3) Please contact our customer service when you cannot receive free-repair service at listed distributor due to the product being given as a gift.
  - 4) This warranty does not apply:
    - to damage caused by misuse, repair, alteration,
    - to damage caused by relocation of storage, transport, drop after purchasing
    - to damage caused by liquid contact, fire, earthquake or other external cause, to incorrect use with abnormal voltage or an undesignated power source (voltage, frequency)
    - shipping damage or damage caused by mishandling by the freight carrier
    - to damage caused by use with other than home use (Example: business purpose)
    - not provide this warranty
    - not provide proof-purchase day, distributor name, your name and reformatted in this warranty
  - 5) Exchange of expendable supplies will be paid.
  - 6) This warranty is for only in Japan.
  - 7) We do not reissue this warranty.
- Purpose for usage of Personal Information and other
- ※ We will maintain and use customer information (copy of warranty) in accordance with limited warranty of free-repair and safety inspection.
  - ※ This warranty is promise of indication of limited warranty for free-repair. It does not set a limit to the customer's legal rights and obligations for a supervisor assigned by this warranty or other business operators. Please contact our customer service when you have questions regarding repair after the end of the warranty period.
  - ※ Please see the manual for "warranty & after-purchase servicing" after the end of the warranty period.

Haslux Co.,Ltd. Beauty Business Division


1-20-5 Ginza Chuo-ku Tokyo ZIP:104-0061 TEL:03-3562-6116 FAX:03-3562-6117

<http://www.exideal.jp>

Warranty

## Manual / Warranty

# Exidéal

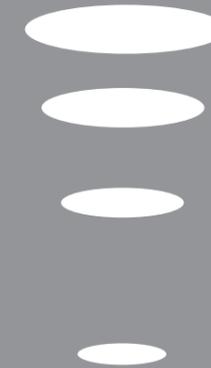

Cover Page	
P. 1-P. 2	Safety and Notes
P. 3-P. 4	Contents of Packages
P. 5-P. 6	Exideal Features
P. 7-P. 8	Exideal Operation Method
P. 9-P.10	10 care support / treatment
P.11-P.12	Wake Up Mode / Q & A
P.13-P.14	Troubleshooting After-purchase servicing
Back Cover Page	Exideal Warranty

Made in JAPAN



Thank you for purchasing our product "Exideal."  
Please read this manual carefully and use our product.  
The warranty is included on the back side of this manual. Please keep it safe.

Haslux Co.,Ltd.

## Please read this instruction manual carefully and use the product safely and correctly.

- These notes provide instruction on the safe and correct use of our product and to prevent injury and damage.
- Below are the classifications for these notes;

### Figures & Classification

-  **Danger** This symbol Indicates an imminently hazardous situation which, if not avoided, may result in death or serious injury.
-  **Warning** This symbol indicates that the user risks serious injury or death if the warning is ignored and the product is operated incorrectly.
-  **Caution** This symbol indicates a potentially hazardous situation, which, if not avoided, may result in minor or moderate personal injury and/or damage to the equipment.

### Explanation of Figure

-  Do not disassemble  This symbol is displayed for acts that require attention or are prohibited. The left symbol illustrates a particular prohibited action such as do not disassemble, repair or modify it, and "Do not disassemble".
-  Prohibit  This symbol Indicates a mandatory action that must be followed. The left symbol illustration "Prohibit" (must be followed explicitly).

 <b>Danger</b>	
	Persons who use a heart pacemaker or electro-medical apparatus that may cause calibration changes should not use this device.

 <b>Warning</b>	
	Please use the exclusive goggles and close your eyes during irradiation. Please do not use anything other than Exideal exclusive goggles.
	<p>Please do not use this product on the face or body for any purpose other than beauty treatment, as other use may damage this product.</p> <p>Please avoid using in the following which, if not avoided, could result in serious injury or skin problems.</p> <ul style="list-style-type: none"> <li>●Epileptic Seizures from the flashing light ●Cardiovascular Disease ●Acute Disease</li> <li>●Allergen ●Boundedness ●Sensitive Skin ●Hypersensitivity ●Abnormal Blood Pressure</li> <li>●Hemophilia ●Infectious Diseases ●Tuberculous Disease ●Under the influence of alcohol or other drugs ●Excessive sunburn and inflamed skin ●Synthetic fibres for sensitive skin ●During the period ●Pregnancy ●Severe fatigue</li> </ul> <p>Please do not use on the following body parts</p> <ul style="list-style-type: none"> <li>●Pimples, Suppuration ●Other than the specified site ●Scars ●Itchiness, burning sensation or freckles from skin disease ●Facial melasma ●Dermatitis caused by skin care product ● Mucous membranes in the mouth ●Excessive sunburn</li> </ul> <p>Please do not use a damaged code, adaptor or power supply. It could cause a fire, short circuit, or deliver an electric shock.</p> <p>Please do not damage, bend, pull, stretch, kink, or place heavy objects on power cords. It could cause a fire, short circuit, or deliver an electric shock.</p> <p>Please use electricity rated at 100-240V power code. Also, please use the Exideal dedicated AC adaptor. Not doing so could cause a fire, short circuit, or other damage.</p>

 <b>Warning</b>	
	<p>Please remove glasses and contact lenses. Wearing them could cause an accident or unexpected problem.</p> <p>Please remove earrings, necklaces, and any kinds of accessories. Wearing them could cause an accident or unexpected problem.</p> <p>Please do not place this product in the bathroom, near too much moisture, or next to water (lavatory vanity). It could be a fire hazard, short circuit, or cause insulation deterioration.</p>
	<p>When you stop using the device please turn off the power and pull the power cord from the outlet.</p> <ul style="list-style-type: none"> <li>●Smoke from this product ●Unusual heat ●Unusual noise ●Unusual smell, etc.</li> </ul>
	Please do not modify or re-assemble without the assistance of a repair engineer. It could result in injury or fire.
	Please do not expose to water. It could cause a short circuit or result in an electric shock.
	Please do not remove AC adaptor with wet hands. It could cause a short circuit or result in an electric shock.
	<p>This device contains extremely high infrared irradiation pulse programs. Operation of TV remote controls or AC is not recommended while using the infrared.</p> <p>Please keep out of the reach of children. It could cause an accident or other problems.</p> <p>Please set the operating time for a short duration if you are a first time user or have sensitive skin. Also check your skin during the session to avoid an accident or other possible problems.</p> <p>Please follow the directions for usage time and frequency of use. Not doing so could cause an accident or other problems.</p> <p>Please stop using if you feel any pain or illness. It could cause an accident or other problems.</p> <p>Please be sure to do a patch-test before using special toner.</p>

 <b>Caution</b>	
	Please avoid using this product the same day you do hair removal or shave the targeted parts of the body. It could cause swelling or redness.
	<p>Please hold the base of the unit in your hands when the LED irradiation panel is angled over 90 degrees. Otherwise, it is dangerous and could fall on your face. Also, please do not try to use this device higher than your face as it is very dangerous to open the LED irradiation panel more than a 90 degree angle.</p> <p>Please do not damage this product or the AC adaptor. It could cause an accident or other problems.</p> <p>Please do not place items on this product or panel.</p> <p>Please do not pull the cable when removing the AC Adaptor.</p> <p>Please do not place the device in the following locations.</p> <ul style="list-style-type: none"> <li>●Exposed to sunlight ●Wet environments ●Dusty environments ●Ignite environments</li> <li>●Rain and wind environments ●High temperature and below zero temperature environments ●Area that generates chemical and gas ●High-vibration environment</li> </ul> <p>Please be sure to do a patch-test before using special toner.</p> <p>Please see page 4 "introduction of patch-test. Please stop using this device and see your doctor if you have any skin problems due to the patch-test.</p>

※Please keep this warranty in an accessible place for easy reference.  
 ※Misuse and incorrect storage of this device are not covered by the warranty, and repair is subject to charge even during the warranty period.

# Contents of Packages / Names and Functions of Display Parts


## 1 Package Contents This package includes the following contents. Before use, be sure to check contents.




### Contents of Goods

- 1 Exideal Product
- 2 Exideal Dedicated Container
- 3 Exideal Dedicated Goggles (Including Cleaning Cloth)
- 4 Exideal Dedicated Goggles Case
- 5 Exideal Dedicated AC Adaptor

## 2 Names of Display Parts

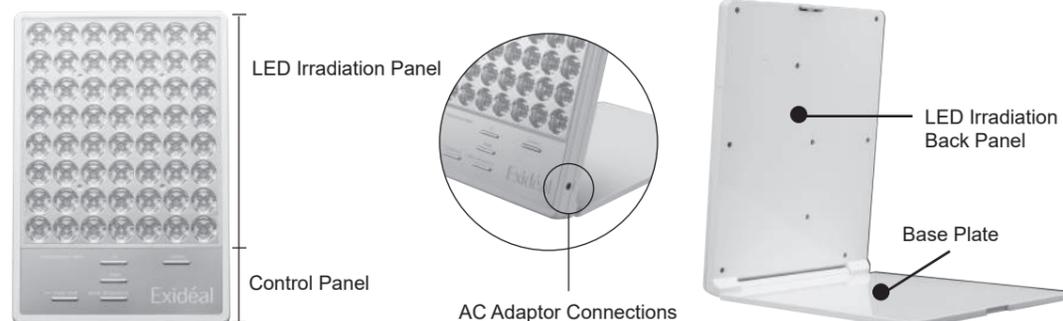


### 【LED Irradiation Panel Names】

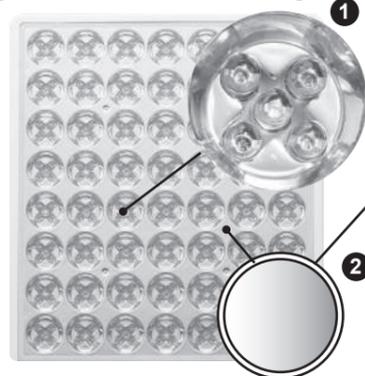

- 1 LED Circle (Buit-in 5 LED)
- 2 Mirror Panel (Mirror surface)

### 【Control Panel Names】

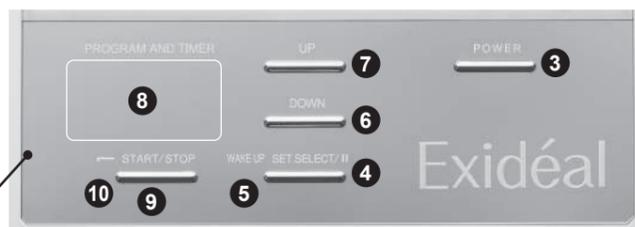

- 3 Power Button (Startup)
- 4 SELECT Button
- 5 Wake Up Setting (Alarm Clock) / Menu Setting / Irradiation Pause
- 6 DOWN Button: Program Selection / Timer Setting
- 7 UP Button: Program Selection / Timer Setting
- 8 Display Monitor (Program, Timer Monitor)
- 9 START Button: Start and Stop Button
- 10 START Light

## 3 Control Panel Function

- 3 [POWER] Button  
Pressing the POWER button (startup) starts all operations with an electronic sound. Pressing the POWER button also ends any ongoing operations.
- 4 [WAKE UP SET SELECT/II] Button (SELECT Button)  
SELECT button has multiple functions.  
(1) Irradiation program and timer selection function  
(2) WAKE UP selection with switching function  
(3) WAKE UP timer selection of hour/minute with switching function  
(4) Program Pause and restart function
- 6 [UP] [DOWN] Button  
UP button is for setting the irradiation program, timer and WAKE UP.
- 7 (1) Pressing the UP button enables you to change to another program number and switch to other programs while choosing an irradiation program. Pressing the UP button once proceeds to the next program. DOWN button returns to the previous program when pressed once. (P01 – P10)  
(2) Press the UP button to change to another control value and to extend the irradiation program timer one minute longer. Press the DOWN button to change to another control value and to reduce the irradiation program timer by one minute shorter. (between 1min and 30min)  
(3) WAKE UP timer can also be used to select hour/minute.
- 8 [Display Monitor]  
(1) Irradiation program: It displays program number on the left "P" and double digits on the right. (Ex: [P:05])  
(2) Timer: It displays double digits of timer "minute" on the right. (Ex: [ 15] "15 min")  
(3) WAKE UP: It displays one digit for minute on the left and double digits for minute on the right (Maximum 9 hours 59 minutes). (Ex: [3:15] "3 hours 15 minutes")  
(4) Ongoing irradiation program: It alternately displays program number and time remaining (by 1 min).  
(5) It displays countdown in 5 second increments for 1 minute when the irradiation program is started.
- 9 [START/STOP] Button (STARAT Button)  
(1) Press button to start program after selecting irradiation program.  
(2) Press button to stop program after starting irradiation program.  
(3) Press button to finish setting WAKE UP after setting WAKE UP time.



## 4 AC电源适配器



- AC: 100-240V
- 50/60Hz
- 24V 1A
- 105.4×44.6×29.2mm

## 5 Dedicated Goggles



※ Goggle use depends on the program. Please read the operation method (P.7) in this manual carefully. Please use goggles when needed.



Dedicated Goggles Case

## Exideal Features

**Exideal uses the 280 LED for high brightness light and LED output driver 15W (watt).4 wavelength high power and Large-Face-One Design beauty device.**

### 1 Easy Use (Only use POWER button and START button for everyday use)

It chooses a treatment program for specific use. There are not so many buttons (4 buttons) and it is easy for everyone to use. The sound effects help guide you.

### 2 Contactless Beauty Device


Treatment is possible without contacting the skin. Contactless treatment dramatically lowers stress on the skin. Also, irradiation does not transfer bacteria like touching.



Maximum 10000 times / second pulse output  
2 times more than regular direct-current instant brightness irradiation

### 3 LED for High Brightness Light

The 280 of 4 wavelength LED for high brightness light is designed for wide-range irradiation home use.

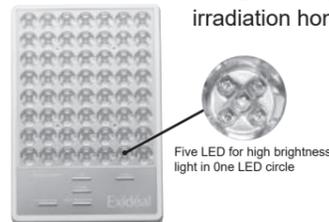

LED circle 56 x 5 LED = **280 LED**  
for high brightness  
**for high brightness light**

950nm wavelength LED (Infrared light): 28  
640nm wavelength LED (Red LED): 112  
590nm wavelength LED (Yellow LED): 112  
460nm wavelength LED (Blue LED): 28

### 4 Maximum LED Output Driver 15W


Outstanding LED output driver for home use.

### 5 Our original optical design

The entire base panel is our original optical design with a finished mirror surface of dermis interflexion. This optical design makes photo-treatment more effective and increases the irradiation rate 120% and 130%.



Base entire panel is our original optical design with a finished mirror surface

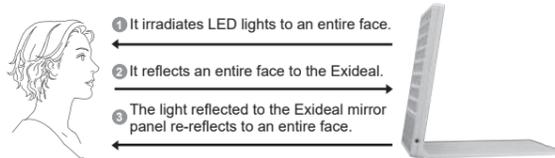

### 6 No stress with hands-free device

You can relax and take care of your skin with this hands-free device. Also, more effectively massages your face during irradiation.



You can massage your skin while receiving a treatment.




You can talk on the phone during a treatment.

#### [Caution]

Please have a base spray in your hands when the LED irradiation panel is angled over 90 degrees. Also, please do not try to use this device higher than your face. It is very dangerous to open the LED irradiation panel more than a 90 degree angle.

### 7 4 wavelength high power and Large-Face-One Design beauty device.

It is the World's First large size device to irradiate your entire face. (Just face size) It covers your entire face and is a very efficient and effective treatment. It possesses a large size design so you can perform a photo-treatment on all of your body parts. This beauty device is hands-free and easily treats out of reach your back. Just Face Size also contains "T-zone treatment."

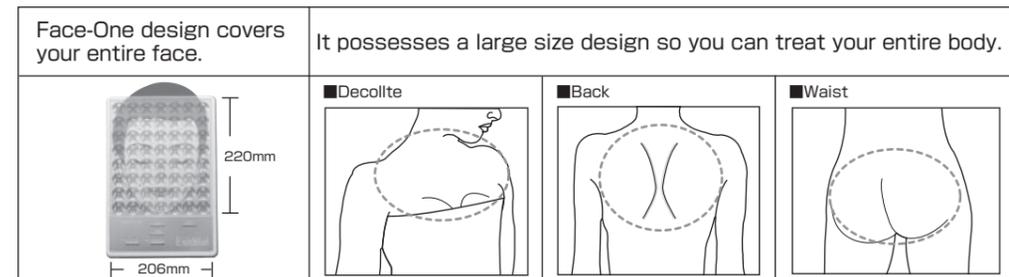

### 8 "Oily T-Zone" Care—Support-Program (P08)




T-Zone is independently controlled so you can treat your entire body and provide T-Zone care at the same time.

• "T-zone" care-support-program's number is "P08."

### 9 MY Memory Function

The machine system memorizes the most recent program (function) used. It displays the previous program (timer) for the next session. You can also make your own original program. Also, the MY memory function memorizes the program forever, even if you remove the adaptor.\* Maximum 30 minutes

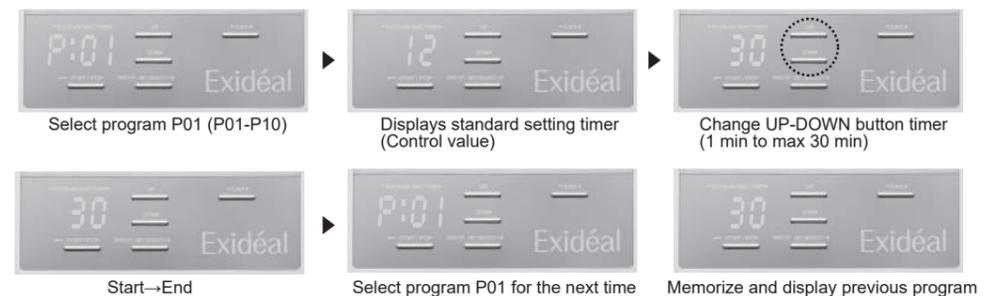

### 10 Wake Up Mode

This function irradiates close to natural light. The timer can be set to use high-drive-output and wide-range irradiation to help you recognize morning light. The morning lights helps to establish routines (like going to bed and waking up early).

#### [Caution]

Please have a base plate in your hands when LED irradiation panel is angled over 90 degree. Also, please do not try to use this device higher than your face and it is very dangerous to open LED irradiation panel more than 90 degree angle.



## Operating Procedures

**1**



**Clean your face with the cleanser and spread the lotion over your entire face.**

- You can use the lotion before, during, and after irradiation.
- \* You can use basic skin-care products such as lotion and cream.

**2**



**Please connect the AC adaptor into the device.**

**3**



**Press a button on the control panel to turn on power**

- (1) Go off a startup-sound.
- (2) Program number on display monitor. (First time sound is P01. Last operation of program number is shown after the 2nd time.)

※ No sound during silent mode.

**4**



**Please select the program you like by checking the program number with the up-down button from the control panel.**

- (1) Select irradiation program with UP/Down button.
- (2) P button: One click and show P01 to P10 → P01 in order.
- (3) DOWN button: One click and show P01 to P10 → P09 to P01 in order.
- (4) Displayed-Program without select program starts immediately when you press the start button.

**5**



**Press SET.SELECT button for time setting ⑥ when you decide to use the program.**


**6**



**Please set the timer (maximum for 30 mins) with the UP/DOWN button from the control panel while checking displayed-monitor for irradiation time.**

- (1) The initial display shows for default value program (Standard irradiation time). It can be operated by procedure ⑦ if there is no change.
- (2) The UP/DOWN button can change the irradiation time.
- (3) UP Button: One click and show 1 to 30 → 1 in order in every 1 mins.
- (4) DOWN Button: One click and show 1 to 30 → 29 to 1 in order in every 1 mins.
- (5) It memorizes previous action time when you press the start button during set to irradiation time.
- (6) It changes configuration of previous action of irradiation program time during irradiation time.
- (7) Program number and selecting time on the display-monitor can be confirmed. To do so press the set-button with set to program and time in a series of operations without the UP/DOWN button.

**7**

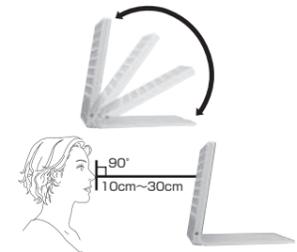

**Adjust for irradiation angle and irradiation distance.**

**8**



**Please wear the goggles if the other programs require the goggles.**

Please use Exideal exclusive goggles when you use this product.

- Programs that require the goggles P01, P03, P05, P07, P08, P09, P10


**9**



**Start irradiation by pressing the start button. Please press the start button once again to abort the program.**

- (1) After you finish setting the irradiation program and timer start the program by using the start button. Pressing the start button will activate the lights on the start-display (red-color).
- (2) Irradiation starts after the 5 second countdown on the display-monitor when choosing a program that uses goggles. (Please place on the goggles within the 5 seconds)
- (3) Pressing the button starts irradiation for programs that do not require goggles.
- (4) It shows on display monitor in 5 seconds during program.
- (5) Pressing the SELECT button during a program stops irradiation with a pause sound. It confirms pause state and changes the start-display to a flashing (red-light).
- (6) Pressing the SELECT program once again restarts the program.
- (7) Pressing the start button during program and pause stops the program in process. It becomes possible to reset the program and also resets the countdown timer.
- (8) Pressing the start button memorizes the executed program. It always memorizes the previous program and shows it in the next selection program.

**2**

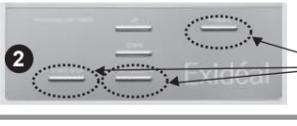

**Use Silent Mode**

- (1) Press POWER button, SET, SELECT / II button, START / STOP button at same time.
- (2) Press START/STOP button to begin irradiation.

※ Please do (1) to stop silent mode.

## 10 Care Support Programs

### Basic irradiation program

Our original skin care support program is effectively combined with the 4 wavelength LED. Please choose your skin type program.

### Caution

Please do not use it if directed by a physician or while taking medication. Our care-support-treatment program is not a medical treatment for the skin.

## Program 01 & 02: Texture firmness skin care support program



**Please use for texture firmness skin care support.**

### Advice and Tips

It is important to moisturize your skin for texture firmness skin care. It is better to do our treatment program after warming your skin with a hot towel and conditioning your skin with a toner, lotion and cream. It treats all age types for texture firmness aging care.



### Using the goggles

For the program P01, please use the exclusive goggles and close your eyes during irradiation. For the program P02, the LED light is not strong, but please close your eyes during irradiation. Program 02 has a time division control, but it has the same power and effect as program P01.

### [P01 With Goggles]

### [P02 Without Goggles]

※ Please close your eyes when using

	Control value Basic irradiation time <b>5 minutes</b>		Control value Basic irradiation time <b>12 minutes</b>
---	--	---	---

## Program 03, 04: Adolescent Skin Care Support Program



**Please use it for smooth skin care support.**

### Advice and Tips

Stand still for aging skin. Please clean your hands when you clean your skin. Please use our program after conditioning your skin with a toner. It is important that you not only take care of your skin, but also take care of the inside of your body by living a healthy daily lifestyle (diet and enough sleep).



### Using the goggles

For the program P03, please use the exclusive goggles and close your eyes during irradiation. For the program P04, the LED light is not strong, but please close your eyes during irradiation. Program P02 has a time division control, but it has the same power and effect as program P03. Caution: Please do not use it vigorously when you have inflammation of the skin.

### [P03 With Goggles]

### [P04 Without Goggles]

※ Please close your eyes when using

	Control value Basic irradiation time <b>5 minutes</b>		Control value Basic irradiation time <b>10 minutes</b>
---	--	---	---

### Caution

Please do not use it vigorously when you have inflammation of the skin. Also do not use it if directed by a physician or while taking medication.

## 10 Care support Programs

### Program 05, 06: Aging and dry skin care program

P05

**Our product prevents dry skin which can be caused by AC or the dry air of winter time.**

**■Advice and Tips**

"Change living environment" and "balance of oil and sweat" are closely related to each other. Dry air and AC cause dry skin throughout one year. Slightly massage skin with oil and cream for counteract skin after conditioning your skin with a toner and serum.

P06

**■Using the goggles**

For the program P05, please use the exclusive goggles and close your eyes during irradiation. Program P02 has a time division control, but it has the same power and effect as program P05.

[P05 With Goggles]		[P06 Without Goggles]	
※Please close your eyes when using			
	Control value Basic irradiation time <b>5</b> minutes		Control value Basic irradiation time <b>10</b> minutes

### Program 07: Beautiful skin care support program

P07

**Please use our product to treat your skin for overexposure to ultraviolet rays.**

**■Advice and Tips**

Your face receives the most exposure to ultraviolet rays and this is a cruel environment. Applying different sunscreens for the different seasons is very important. Use the beauty ingredients contained in the toner and serum everyday. Be sure to practice a healthy daily lifestyle (proper diet and enough sleep).

[P05 With Goggles]	
※Please close your eyes when using	
	Control value Basic irradiation time <b>5</b> minutes

※The LED irradiation is extremely high, please use the exclusive goggles.

### Program 08: The T-zone program helps to treat oily skin.

P08

**Please use our product for problems like "Annoyed T-Zone"**

**■Advice and Tips**

T-zone becomes oily and hard skin when unbalanced between sunburn and fluid. It is caused by oxidation when skin suffers excessive sunburn which prevents water evaporation dryness and fluid balance. Be sure to clean sunburn to prevent skin problems.

[P08 With Goggles]	
※Please close your eyes when using	
	T-zone
	Control value Basic irradiation time <b>5</b> minutes

### Program 09: Total care support in the fall and winter

P09

**Basically, total care support for the fall and winter time. Multiple functions of support care in short programs. The care support program is a mix of the short programs P01, P05 and P07.**

**■The contents of multiple functions**

- (1) P01: Skin texture and care support.
- (2) P05: Aging and dry skin care.
- (3) P07: Healthy skin and a smile! Beautiful skin care support.

※It conducts all three programs automatically every twenty seconds.

**[P09 With Goggles]**

※Please close your eyes when using

	Control value Basic irradiation time <b>8</b> minutes
---	---

### Program 10: Total care support in the spring and summer

P10

**Basically, total care support for the spring and the summer time. Multiple functions of support care in short programs. The care support program is a mix of the short programs P01, P07 and P08.**

**■The contents of multiple functions**

- (1) P01: Skin texture and care support.
- (2) P07: Healthy skin and a smile! Beautiful skin care support.
- (3) P08: Helps to treat the oily T-zone.

※ It conducts all three programs automatically every twenty seconds.

**[P10 With Goggles]**

※Please close your eyes when using

	Control value Basic irradiation time <b>8</b> minutes
---	---

### Cautions on Use / Safety Functions

**■Cautions on Use**

- 1 Please use this device at a unit temperature between -15°C and +40°C. Adequate irradiation treatment might not be received if the temperature is not within normal range of use. It could require the operation safety function of the micro computer. We recommend that you use it below 30°C. Please do not use in a sauna. It may cause deformation and failure.
- 2 Please do not use the device in difficult environments like wrapping the device in a blanket, as the device cannot radiate. Adequate treatment may not be received if the operating safety function is activated near a stove or AC outlet where high temperatures can result in burns or deformation of the unit.
- 3 Please avoid using this device in places of condensation such as a bathroom or washbasin which can cause malfunction or failure of the micro computer. Also, please do not use it in the bathtub as it can cause electric shock.

**■Safety Functions**

- 1 This device has monitoring functions in which the micro computer can detect the internal temperature and lower output automatically. Please be cautious and use carefully in appropriate environments. Excessive use, a very hot environment, or radiation distraction can cause irradiation to stop operating and error codes "E-01" or "E-3" to be displayed. Please press the POWER button for restarting. Please remove AC adaptor and POWER OFF if it still displays "E01" after restarting. Then please try again after one hour.
- 2 This device has monitoring functions which detect the internal temperature and it also has safety equipment. It stops automatically when the device temperature rises abnormally from a stove or other heat source. Please remove AC adaptor and POWER OFF for one hour when it displays "E02." Also avoid exposing the device to abnormal temperatures when using it.


**A lack of sleep and unstable lifestyle are enemies of beautiful skin.  
“Beauty and Healthy” is a theme.**

**We support a regular healthy routine and lifestyle with light in the morning.**

Our original lights help to establish healthy routines. Automatically turn our original highlights on by setting the timer and it will alert you in the morning. The lights help to establish healthy routines like going to bed and waking up early.


### [Functional Overview]

The built-in timer operates for twenty-four hours. For example, if it was 11 p.m. now and the timer was set for 7 a.m., it means the timer goes off eight hours later. After eight hours, it automatically begins irradiating which makes your morning refreshing. The function iterates everyday as long as you do not release the device or unplug the device.

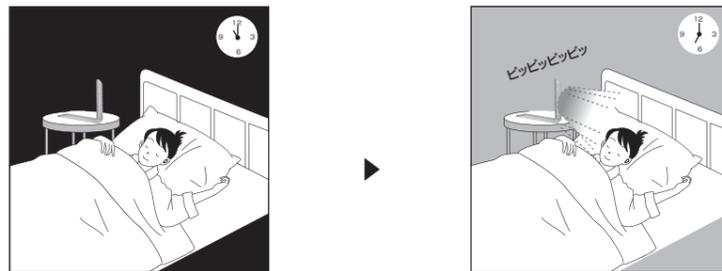

### [Wake Up Mode Operating Procedures]

- (1) Press POWER button to turn on power.
- (2) Press SELECT button for 2 seconds when irradiation program is not ongoing. It flashes WAKE UP next to the SELECT button when “9:00 (control value / after 9 hours)” is on the display monitor. Also it flashes “9 hours” on the left at same time. Please set time with the UP/ DOWN button (0-9).
- (3) Press SELECT button again and it flashes “0 (10 min)”. Please set by 10 minutes with UP / DOWN button (0-9).
- (4) Press the SELECT button again and it flashes “0 (1 min)”. Please set by 1 minute increments with the UP / DOWN button (0-9).
- (5) Please press the START button when you finish setting (3) to (5). A “beep” sound will be heard once the wake-up program is set. It flashes WAKE UP next to the SELECT button on the display monitor.
- (6) The same program will operate every day at the same time once wake-up mode is set. If you would like to cancel the program, remove the AC adaptor or cancel wake-up mode. It flashes WAKE UP on the display monitor when wake-up mode is set. Turn off the POWER button and wake-up mode flashes on the display monitor.

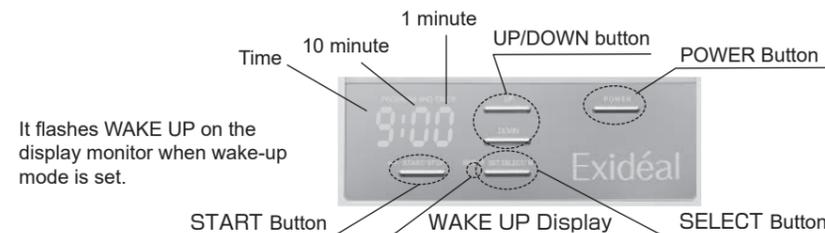

### [WAKE UP Mode Operation]

- (1) Wake-Up mode startup begins the wake-up irradiation program for 15 minutes.
- (2) Pressing the START, SET or SELECT button will POWER OFF irradiation during the wake-up mode program.
- (3) The unit will POWER OFF automatically after it finishes irradiating for 15 minutes.

### [WAKE UP Mode Termination]

- (1) Press POWER button to turn on power.
- (2) Press SELECT button for 2 seconds. (return to select mode)
- (3) Press the SELECT button again for terminating the program and it will flash WAKE UP on the display monitor.  
※A power outage will terminate wake-up mode. If a power outage occurs, please reset the device.  
(Remove AC adaptor for termination)

### Q&A

Frequently Asked Questions	
Q1	<b>Q:</b> What is the LED? How do the lights effect the skin? <b>A:</b> LED is an abbreviation for Light Emitting Diode. A color (wavelength) is characteristic of the LED light which is a semiconductor used in home appliances and tags. High value wavelength strengths have been carefully chosen for skin care. Each LED has its own color that is created by a specific periodic wave called a wavelength. This wavelength is very important for skin care.
Q2	<b>Q:</b> Does it affect my eyes? <b>A:</b> NO. It does not directly affect your eyes. There are stronger lights in our daily life. Some programs emit a momentary brightness that can blind you temporarily, please use the goggles required by those programs. (*1) Please close your eyes even if the other programs do not require the goggles.
Q3	<b>Q:</b> Is it possible to irradiate an entire body? <b>A:</b> No problem..Please use the goggles required by those programs(*1) and also do not look directly at the irradiation. Exideal is a large-face-one design so you can use it for an entire body.
Q4	<b>Q:</b> I have a problem with acne. Is there a bad program for acne? <b>A:</b> Please try the adolescent program. Please skip treating your skin if you experience bad irritation.
Q5	<b>Q:</b> Which programs are best for daily use? Is it a problem if I use two or three different programs? <b>A:</b> You do not have to switch the programs because it switches them automatically when you use total care support for the spring and summer, or the fall and winter. You can switch the programs on your own. The ultraviolet is not harmful so enjoy using it. Please do not treat the same area for more than 30 minutes.
Q6	<b>Q:</b> Is there any irritation during a session? <b>A:</b> Is it a problem if I use two or three different programs? This device is contactless and does not pass a current. This eliminates electrical stimulation and makes the device very safe. It can cause skin to dry a little so please moisturize your skin.
Q7	<b>Q:</b> Is it okay to use the device after my skin is irritated from tanning? <b>A:</b> Please do not use it for three days when your epidermal skin is burned by ultraviolet light. You can use it after cleaning your face as long as there is no stiffness or any other trouble with your skin.
Q8	<b>Q:</b> Is there any difference in effectiveness between a short distance and a long distance? <b>A:</b> Effectiveness and results can differ if the lights are weakened by being far from the LED irradiation. Please use it at distance between 10cm and 30cm.
Q9	<b>Q:</b> Is wake-up-mode a fire hazard if I forget to turn it off during a trip? <b>A:</b> If you are out of town or do not use the device for a long period of time, please unplug the AC adapter.
Q10	<b>Q:</b> How long can I irradiate a light on my face? <b>A:</b> Irradiating for a long period of time is no problem. But it takes time to make your skin beautiful. Our device helps you to make your skin beautiful so you can keep treating your skin for daily use.
Q11	<b>Q:</b> How much does using wake-up-mode for many months potentially cost me? <b>A:</b> Using wake-up-mode twice a day can cost you approximately 39 yen a month.
Q12	<b>Q:</b> How do I care for the device after I finish treating my skin? <b>A:</b> The device is contactless and requires minimal care. Please use detergent and soft cloths to clean the irradiation surface area and surrounding buttons. Please do not use alcohol, thinner or nail polish remover for cleaning as they can damage the device.
Q13	<b>Q:</b> Can I purchase spare parts if needed? Can the device be repaired? <b>A:</b> You can do both. Go to the shopping cart. Please contact the shop where you purchased the device or our Exideal customer service.
Q14	<b>Q:</b> Can I use it in other countries? <b>A:</b> Yes. Exideal is compatible with 100-240V, please check the voltage of other countries. If you can change the shape of an outlet, you do not need to have a voltage converter.

## Troubleshooting

Trouble Condition	Reasons	Remedy
<b>“E-01” error message on display</b>	This error can result from excessive use, a very hot environment, or because of a radiation distraction	Please press POWER button for restarting. Please remove AC adaptor and POWER OFF if it still displays “E-01” after restarting.
<b>“E-02” error message on display</b>	The device stops automatically when its temperature rises to an abnormal level from a heat source like a stove.	Please remove AC adaptor and POWER OFF for one hour when it displays “E-02.” Also avoid using the device in abnormal temperature conditions.
<b>“E-03” error message on display</b>	This error can result from excessive use, a very hot environment, or because of a radiation distraction.	Please press POWER button for restarting. Please remove AC adaptor and POWER OFF if it still displays “E-03” after restarting.
<b>Turn on but nothing displays.</b>	The AC adaptor is not connected properly into the device.	Please connect AC adaptor for sure.
<b>Wake-up mode is disabled.</b>	The AC adaptor is not connected properly to the device or there was a power outage.	Please connect AC adaptor for sure and reset wake-up mode for power outage.
<b>The LED irradiation does not operate and the display disappears.</b>	Was there any static discharge by a person during the winter?	This device has micro computer monitoring functions that detect large static discharges. When detected, the device resets itself. Please remove the AC adaptor and POWER OFF for five seconds. Then restart it.
<b>It does not operate with the TV, AC or remote control.</b>	There might be an infrared influence between the device and the TV receiver, AC or remote control.	Please change direction of the irradiation surface to protect the light from the TV, AC, or remote control. Turning off those devices would also eliminate the infrared influence.

**[Caution]**

Please contact our customer service if this troubleshooting guide does not answer your questions.

## [After-Purchase Servicing]

**■Warranty**

The warranty is attached to the back cover. Please fill in all necessary information such as location purchased and the date of purchase, and keep it safe. Please read this warranty description carefully.

**■Warranty period**

The guaranteed term is one year from the date of purchase.

**■Request for repair**

Please contact our customer service if you still have problems after troubleshooting. Even if you referred to your user's manual and removed the AC adaptor per the detailed instructions.

※Shipping costs shall be borne by the Customer. We pay shipping costs after repair.

**●Warranty period**

We will repair without charge according to the regulations of warranty. Please bring or send proof-of-purchase. Proof-Of-Purchase is required to verify the warranty period with our customer service.

**●After warranty period**

When the product still performs, we will repair with charge based on the needs of a customer. Please contact the original distributor or our customer service.

### Accessory Price List

AC Adaptor	5,000yen (without tax)
Goggles	3,500yen (without tax)
Original Container	3,000yen (without tax)

## [Specifications] MADE IN JAPAN

**■Product**

Power Supply Type	Accessory Dedicated Adapter	Rated Input Power(Power Source)	DC24V ----- 1A
Power Dissipation	20W (Maximum)	Standby Power Consumption	Under 3W
LED Output Drive	15W (Maximum)	LED Output Wavelength	420~970nm
Weight	1.66Kg	Optimum Temperature	Operating Temperature Range /-5°C~+40°C
Product Dimension	Width:206mm Length:300mm Height:30mm		Storage Temperature Range /-5°C~+45°C

\* This product is a beauty device and is not for medical treatment.

## Haslux Co., Ltd. Beauty Business Division

**●Customer Support**

ZIP:104-0061  
1-20-5 Ginza Chuo-ku Tokyo  
TEL:03-3562-6116 FAX:03-3562-6117

**●Osaka Merchandise Center**

ZIP:566-0055  
2-2-6 Shinzaike Settsu-city Osaka  
TEL:06-6827-0307 FAX:06-6827-0304